

### INTENDED USE

Healeon HD PRP is designed to be used for the safe and rapid preparation of autologous platelet-rich plasma (PRP) from a small sample of peripheral blood at the patient point-of-care. The PRP is mixed with autograft and/or allograft bone prior to application to a bony defect for improving healing characteristics.

### PRODUCT DESCRIPTION

A blood collection tube in a peel pouch or blister tray kit used for Platelet Rich Plasma (PRP) preparation from peripheral blood. The product is provided sterilized and non-pyrogenic. The product contains ACD-A anti-coagulant and is to be used only for the preparation of PRP with Healeon HD PRP.

### NOTE:

- For Single Patient Use Only. Once used, the item must be discarded following the protocol for biohazardous waste disposal, as required by the facility where it was used.
- The PRP needs to be centrifuged within 6 (six) minutes of the blood being drawn into tube.
- The PRP needs to be used within 4 (four) hours after being processed.

### STORAGE

Store tubes at 4-25° C unless otherwise noted on the package label.

### CAUTION

- PRP must be used on the same patient from which blood was drawn.
- Do not use PRP tube if any foreign matter is visible.
- Do not use tubes that are, or appear to be cracked or damaged.
- Do not use PRP tube after its sterilization expiration date.
- Do not re-sterilize or attempt to clean the PRP tubes.
- Once prepared, use PRP within 4 (four) hours after drawing blood.
- Practice Universal Precautions to protect from blood spatter, blood leakage to minimize the potential exposure to blood-borne pathogens.
- Handle all biologic samples and blood collection devices per the procedures established by your facility.
- Appropriate medical attention is required in the event of any exposure to biologic samples.

### REQUIRED ITEMS NOT PROVIDED FOR SPECIMEN ACQUISITION AND PROCESSING

Have these items assembled and ready prior to drawing blood:

- Centrifuge capable of generating 1600 G Force.
- Appropriate items to draw the blood specimen.
- Use gloves, eye protection, gowns, other apparel and disposable items to protect personnel from exposure to blood-borne pathogens or other potential harmful materials. Follow the policies and procedures of your facility for venipuncture and blood draw procedures.
- Alcohol swabs for cleansing draw site.
- Sterile gauze.
- Tourniquet.
- Sharps disposal container for needles, syringes and needle/holder device.
- PRP tube, mixing tube and transfer device.


### GLOSSARY OF SYMBOLS

REF	Catalog Number	LOT	Batch Code	Use By Date
Single Patient Use		Do Not Use If Damaged		Manufacturer
Temperature Limit		Keep Out of Sunlight		Keep Dry
STERILE R	Sterilized using Irradiation	See Instructions for Use		Do Not Restерilize

### PROCEDURE STEPS TO OBTAIN PRP

#### 1) DRAW BLOOD

Remove PRP tube from pouch. Wipe top with clean alcohol swab prior to blood draw. Following accepted blood collection guidelines, draw the patient's blood into PRP tube.

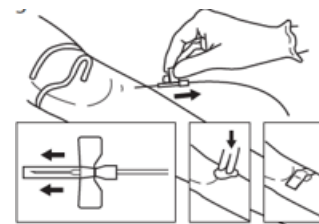

#### 2) INVERT TUBE

Invert PRP tubes 7 (seven) times to mix anti-coagulant with patient blood. Wipe tops with clean alcohol swab. Centrifuge tube within 6 (six) minutes of blood draw.

#### 3) SPIN PRP

Turn on centrifuge and unlock lid. Place PRP tubes into individual tube holders across from each other to achieve a balanced load in the centrifuge, or use appropriate counter balance. Lock the lid.

### IMPORTANT NOTE!

#### HORIZONTAL OR FIXED ANGLE ROTOR -

- 15 ml tube - Spin for 10 (ten) minutes at 1100 G's
- 30 ml tube - Spin for 10 (ten) minutes at 1600 G's

Please refer to your facility's centrifuge manual for proper use and rotor style.



#### 4) PRP EXTRACTION - INVERT TUBE

Remove tubes from tube holders. Using aseptic technique, connect sterile syringe to a long needle and insert long needle through the stopper; aspirate the PRP into the syringe. Once the PRP is transferred into the syringe, the product is ready for use.

Please refer to Table 1 for appropriate PRP volume.

TABLE 1:

BLOOD DRAW VOLUME	PRP VOLUME
15 ml	7.5 ml
30 ml	15 ml

**CAUTION:** Do not insert the needle into the separating gel. If needle comes into contact with the separating gel, discard the product.

### CENTRIFUGATION

**CAUTION:** Follow manufacturers instructions on the use of the centrifuge at your facility.

- Centrifuge carriers and inserts should be of the size specific for the tube used.
- Always allow the centrifuge to come to a complete stop before removing tubes.
- Contra balance tubes will be required to balance the weight of the tubes during centrifugation. When using a contra balance tube, place the balance tube with the same weight as the filled PRP tube, opposite of the PRP tube.
- Once the centrifugation spin cycle has stopped, open lid of centrifuge and examine tubes for any breakage or leaking prior to use.

**CAUTION:** Do not attempt to remove any broken tubes, pieces of glass or other debris with your hands. Use some type of hemostat or forceps to remove items that are damaged. Follow centrifuge manual provided by the manufacturer for disinfection instructions.

**CAUTION:** Federal Law (USA) restricts this device to sale by or on the order of a physician.



**BIMINI HEALTH TECH**  
8400 Bellevue Drive, Ste 125  
Plano, TX 75024, USA  
Phone: 858-348-8050  
www.biminihealthtech.com